



OPTO-EDU (BEIJING) CO., LTD.

F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China
Tel:+8610 88696020 Fax:+8610 88696085

A15.1091

Polarizing Microscope, Semi-APO



A15.1091 Series Microscope



A15.1091-R
Polarizing Microscope,
Semi-APO, Reflect



A15.1091-T
Polarizing Microscope,
Semi-APO, Transmit



A15.1091-TR
Polarizing Microscope,
Semi-APO, Transmit/Reflect

Product Details



Seidentopf Trinocular Viewing Head
0°~35° adjustable viewing head can be operated in a more comfortable position. Minimize muscle tension and discomfort associated with working long hours.



Intermediate Unit

The intermediate tube incorporates a Bertrand lens as standard, enabling both the observation and capture of conoscopic and orthoscopic images. Precisely and quickly interference images adjustment. The high precision slide-type analyzer can be rotated a full 360° with a precision vernier scale and provides unparalleled image contrast.




Reflect Epi Polarizing Illuminator

Up to 6 observation modules can be placed in the rotating disc structure, while providing full extension capacity can be quickly switched. Easy to use and improve work efficiency with just one touch.

Compensator Slot

On nosepiece there is compensator slot. Therefore, various compensators can be used to enhance the advanced quantitative measurement of the signal of the weak birefringent material.

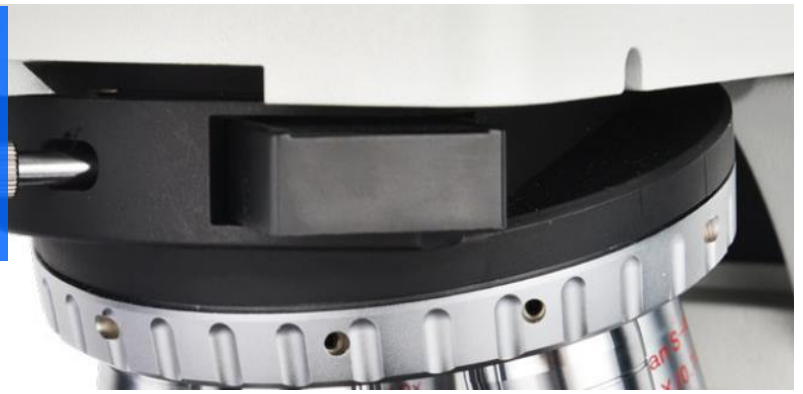

Center Adjustable Quintuple Nosepiece

Up to five objectives can be mounted and all objectives positions are centerable. The nosepiece can be equipped with 6 objectives with different magnification rates to obtain abundant sample information.



A5P.1090 Strain-Free Semi-APO Infinity Plan Objective For Polarizing View

With the deepening of the study of polarized light, the accuracy and repeatability of the data are increasingly required. Objective lens and condenser lens, as important components of the microscope, play a vital role in the imaging effect. As a polarizing microscope, not only does it require high resolution, but also a stress-free optical component, which ensures that the observed birefringence is from the sample and not from the optical element.

High Precision Revolving Round Stage

Diameter 190mm, Center Adjustable, 360° Scale, Minimum Division 1°, Minimum Reading 6' by Means of Vernier Scale, 45° Click Stop Knob. Attachable Moving Stage, Moving Range 30x30mm.



ECO Function

The transmitted light would be off automatically after 30 minutes from operators leave. It can not only save energy, but also keep the lamp life.


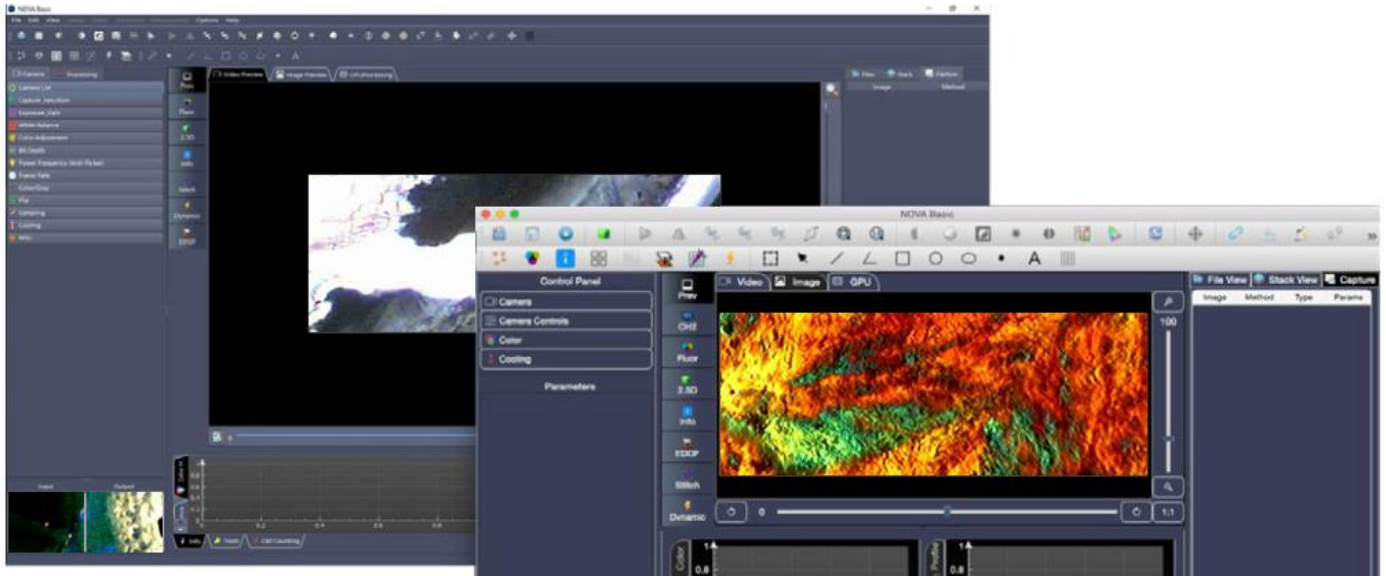


Image Analysis System

Today, the research work environment requires instruments that can be adapted to each person's different workflow. NOVA Basic can also serve as an observation instrument for today's popular operating systems. NOVA Basic microscopic image analysis software allows seamless integration between acquisition, processing, measurement and microscope.

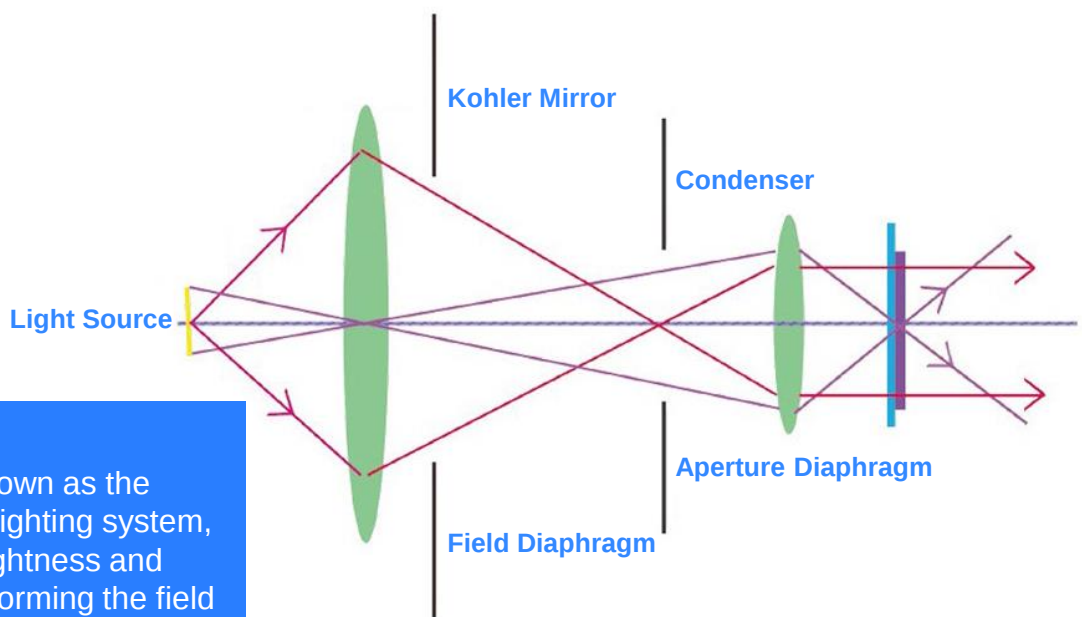


Kohler Lighting


Kohler Lighting is known as the perfect microscopic lighting system, providing perfect brightness and uniformity. While uniforming the field of view, it is possible to expand the functions of the microscope.

Product Application

Traditional polarizing microscopes are used in geology, mineralogy, metallography and geological exploration. However, with the deepening of various disciplines, more and more fields have put forward demands for polarizing microscopes, including not only industrial fields such as construction, material analysis, semiconductor production and detection, but also natural science fields such as forensic science, medical inspection and biological research. At the same time, the application space and extensibility of microscopes are required to be higher, and customized microscopes with higher image clarity, better contrast and higher efficiency are needed. Polarizing microscope NP900 is powerful and comfortable. It is your reliable partner in scientific research and industrial inspection. It can realize many kinds of observation and provide excellent optical imaging quality.

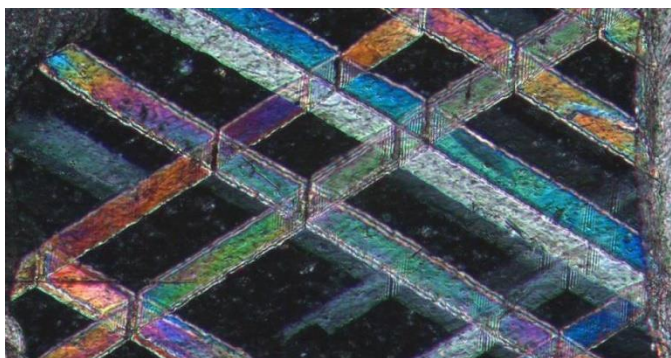

Calcite

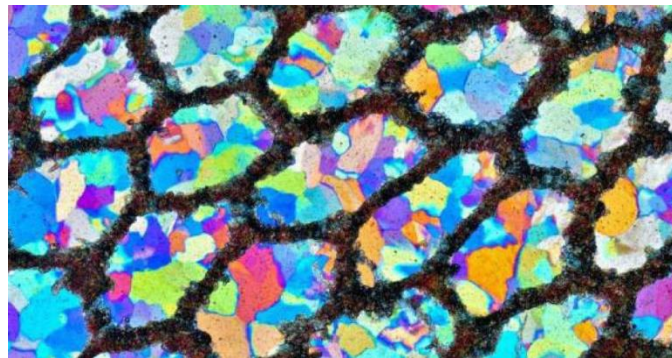

Coral

Accessories



A54.1094-R Moving Ruler

In order to improve the operating efficiency of the microscope, the double-layer mechanical moving ruler can be easily installed on the rotating stage, and the sample can be easily moved. The moving range: 30x30mm

A5P.1090 Multiple Types Of Compensators


A variety of optical compensators are available in A15.1091 that enhance the signal of weak birefringent materials, making them easier to see. At the same time, the compensator can also be used for qualitative measurements to determine multiple horizontal delays.



Software Function



A59.2225

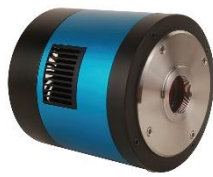


A59.2226

With high speed USB3.0 digital camera and professional image processing software, Opto-Edu microscope can work with computer to get various advanced function done easily. Real Time / Static Measure, 2D Image Scan & Stitching, 3D Depth of Field Fusion, Fluorescent Image Synthesis, Cell Counting and etc


◆ 2D Image Scan & Stitching

By collecting and importing images in real time, the software can quickly stitch together to form a large-size and high-resolution image

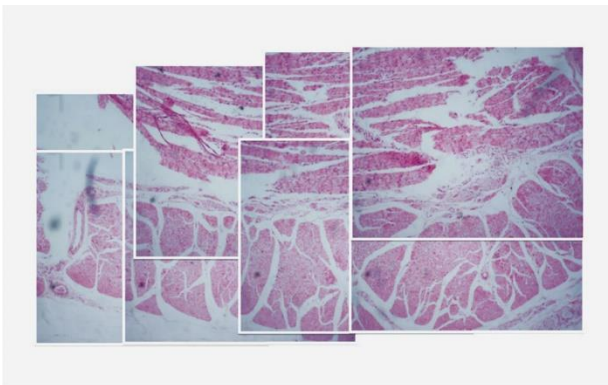

◆ 3D Depth Of Field Fusion

Users can collect multiple images with different focal lengths by fine-tuning the focal length, and synthesize one image for output. Suitable for specimens that require a certain depth of field or poorly made sections

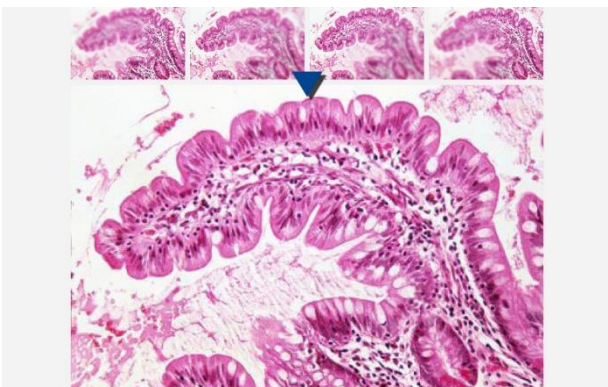


◆ Real Time / Static Measure

In cell observation and section observation, you need to use the measurement function. To determine the cell size, cell gap, synapse length and other data. The software can provide measurement of distance, angle, rectangle, circle, ellipse, etc.




◆ Cell Counting

Customize cell counting requirements, automatically count and count the shape information of cells, including: size, location, volume, circumference, brightness, etc. And all data including processed images can be saved as excel sheet

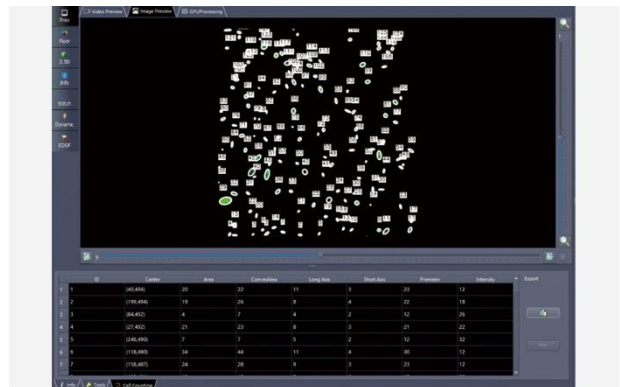

◆ Fluorescent Image Synthesis

By collecting or importing images of different fluorescence channels, users can obtain images after fluorescence synthesis. For the image of each channel, the displacement in the x direction and y direction can be adjusted to achieve the fine-tuning effect

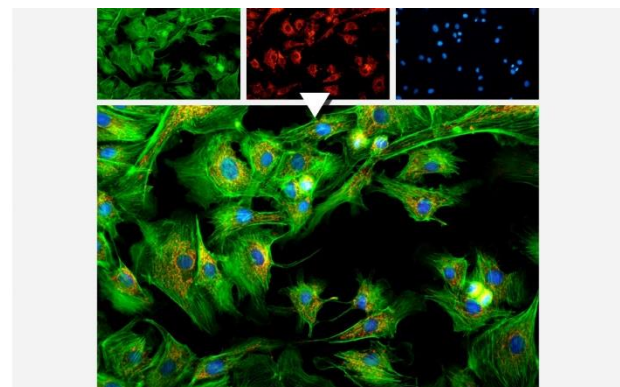


Product Specification


A15.1091 Series Polarizing Microscope		A15.1091			Cata.No.
		-T	-TR	-R	
Optical System	NIS60 Infinite Optical System	●	●	●	
Observation Method	Bright Field	●	●	●	
	Dark Field,	○	○	○	
	Polarizing	●	●	●	
	Phase Contrast	○	○	○	
	Flourescent	○	○	○	
	DIC	○	○	○	
Main Body	BH Manual Body + Halogen Illumination.	●	●	●	A54.1090-BH
	BL Manual Body + LED Illumination.	○	○	○	A54.1090-BL
	ATH Semi-Auto Body + Halogen Illumination.				
	Auto Nosepiece + Auto Condenser + Auto Brightness Adjust	○	○	○	A54.1090-ATH
	ATL Semi-Auto Body + LED Illumination.				
Head	Auto Nosepiece + Auto Condenser + Auto Brightness Adjust	○	○	○	A54.1090-ATL
	Seidentopf Binocular Head, Inclined 30°, Interpupillary Distance 47-78mm	○	○	○	A53.1090-B
	Seidentopf Trinocular Head, Inclined 30°, Interpupillary Distance 47-78mm, 3 Level Ligth Split Switch E100:P0/E20:P80/E0:P100	○	○	○	A53.1090-T
	Ergo Tilting Trinocular Head, Inclination 0°~35°, Interpupillary Distance 47-78mm, 3 Level Ligth Split Switch E100:P0/E20:P80/E0:P100	●	●	●	A53.1090-TT
Eyepiece	SW10x/22mm, High Eyepoint, Diopter Adjustable, Dia.30mm	●	●	●	A51.1090-1022
	SW10x/22mm, With Scale Cross, With Locating Pin, For Polarizing	●	●	●	A51.1090-1022S
	SW10x/22mm, With Cross Hair, With Locating Pin, For Polarizing	○	○	○	A51.1090-1022C
	SW10x/22mm, With Grid Plate, With Locating Pin, For Polarizing	○	○	○	A51.1090-1022G
	EW12.5x/16mm, High Eyepoint	○	○	○	A51.1090-12516
	WF15x/16mm, High Eyepoint, Diopter Adjustable, Dia.30mm	○	○	○	A51.1090-1516
	WF20x/12mm, High Eyepoint, Diopter Adjustable, Dia.30mm	○	○	○	A51.1090-2012
Nosepiece	Manual Nosepiece, Quintuple, Backward, Center Adjustable	●	●	●	A54.1091-5M
	With Slot For Polarizing Compensator Slider Or DIC Slider	●	●	●	
	Protect Cover For Nosepiece Holes	●	●	●	A54.1091-C
NIS60 Strain-Free Semi-APO Infinity Plan Objective	5x/0.15, W.D.20mm, Semi-APO, PL, No Cover Glass	●	●	●	A5P.1091-5
	10x/0.3, W.D.11mm, Semi-APO, PL, No Cover Glass	●	●	●	A5P.1091-10
	20x/0.45, W.D.3.0mm, Semi-APO, PL, No Cover Glass	●	●	●	A5P.1091-20
	50x/0.8, W.D.1.0mm, Semi-APO, PL, No Cover Glass	●	●	●	A5P.1091-50
	100x/0.9, W.D.1.0mm, Semi-APO, PL, No Cover Glass	●	●	●	A5P.1091-100
Working Stage	High-precision Pre-Adjusted Revolving Round Stage, Diameter 190mm, Center Adjustable, 360° Scale, Minium Division 1°, Minium Reading 6´ by Means of Vernier Scale, 45° Click Stop Knob.	●	●	●	A54.1094
	Attachable Moving Stage, Moving Range 30x30mm	●	●	●	A54.1094-R
Condenser	Swing-Out Condenser NA0.9/0.25	●	●	●	A56.1091-S
Focusing	Coaxial Coarse & Fine Focusing, Fine Division 0.001mm, Focusing Range 35mm, Coarse Stroke 37.7mm, Fine Stroke 0.1mm, Can Exchange Hand Wheel Between Left/Righ	●	●	●	

Transmit Light Source	Transmit Kohler Illumination, Brightness Adjustable, 12V100W Halogen, External Lamp House ECO Function Support Auto Power Off After 30 Mins From Operator Leave To Save Energy	●	●	-	A56.1090-12V100W
	Filter Holder On Base, Can Hold 3 Filters	●	●	-	A56.1092-H
	Filter LBD	●	●	-	A56.1092-LBD
	Filter Green	●	●	-	A56.1092-G
	Filter Yellow	●	●	-	A56.1092-Y
	Filter ND6	●	●	-	A56.1092-ND6
	Filter ND25	●	●	-	A56.1092-ND25
Polarizing Attachment	Transmit Polarizer, Attached To Condenser Bottom, 360° Rotatable	●	●	-	A5P.1090-TP
	Reflect Polarizer Slide, Insert In Slot of Reflect Light Source	-	●	●	A5P.1090-RP
	Reflect/Transmit Analyzer Slide, 360° Dial Rotate, Minimum Scale Reading 0.1°	●	●	●	A5P.1090-TRA
	Bertrand Lens Intermediate Unit, With Slot For Compensator Slide	●	●	●	A5P.1090-CB
	λ Compensator Slide	●	●	●	A5P.1090-C λ
	λ 1/4 Compensator Slide	●	●	●	A5P.1090-C1/4 λ
	Quartz Wedge Slide	●	●	●	A5P.1090-CQ
Reflect Polarizing Light Source	Blank Slide	●	●	●	A5P.1090-CS
	Reflect Epi Polarizing Illuminator, Turret Disc With 6 Positions For Filter Block, Kohler Illumination	-	●	●	A5P.1090
	12V100W Halogen Lamp Housing	-	●	●	A5P.1090-100W
	BF1 View Block	-	●	●	A5P.1090-B1
	BF2 View Block	-	●	●	A5P.1090-B2
	Circularly Polarized Light Block	-	●	●	A5P.1090-CP
	Filter Blue	-	●	●	A56.1093-B
	Filter Green	-	●	●	A56.1093-G
	Filter Yellow	-	●	●	A56.1093-Y
	Filter Forested	-	●	●	A56.1093-F
Adapter	Eyepiece Adapter Dia.23.2mm	○	○	○	A55.1090-E
	C-Mount 1.0x	○	○	○	A55.1090-10
	C-Mount 0.5x	○	○	○	A55.1090-05
Other Accessories	Working Stage Holder Bracket	●	●	●	A54.1096
	Adapter To Adjust Eye Position	○	○	○	A54.1096-A1
	Adapter To Lower The Stage Position 1"	○	○	○	A54.1096-A2
	Immersion Oil	○	○	○	A50.1090-01
	Allen Wrench	●	●	●	A50.1090-02
	Power Cord	●	●	●	A50.1090-03
	Short Eye Cover, For Eyepiece	○	○	○	A50.1090-04
	Long Eye Cover, For Eyepiece	○	○	○	A50.1090-05
	Eyepiece Micrometer, Cross	○	○	○	A50.1090-06
	Adapter Ring To Install Eyepiece Micrometer	○	○	○	A50.1090-07
	USB Cable	○	○	○	A50.1090-08


Note: "●" In Table Is Standard Outfits, "○" Is Optional Accessories "-" Is Unavailable

System Diagram & Size(mm)

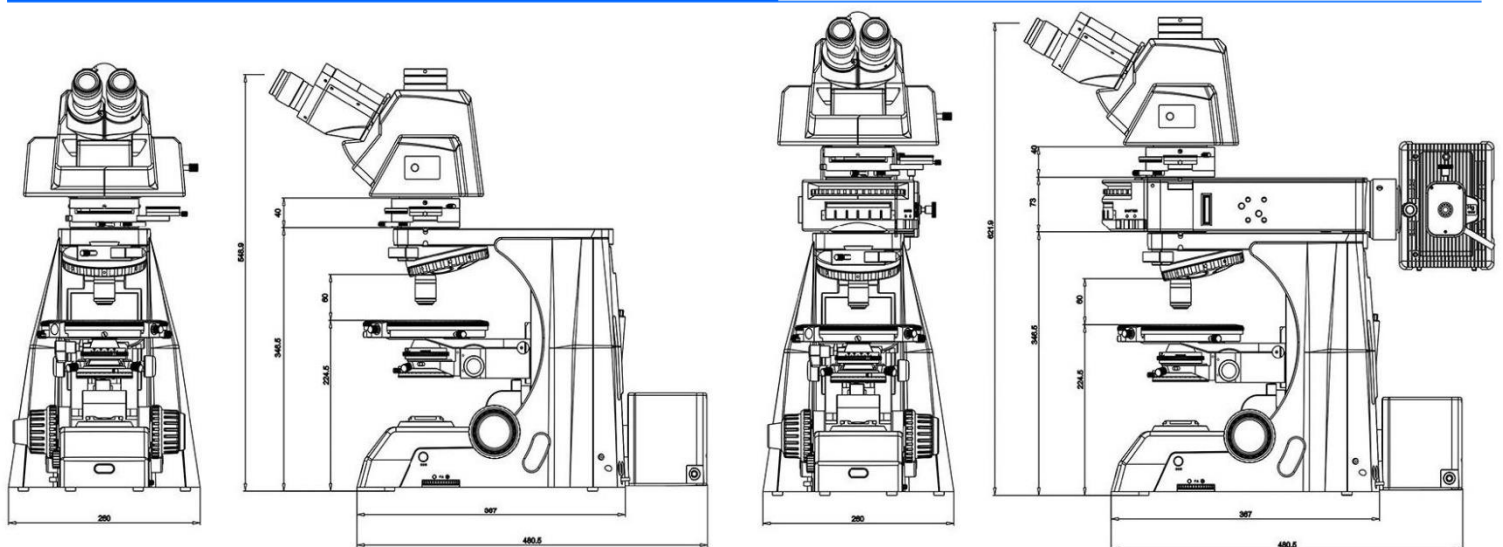


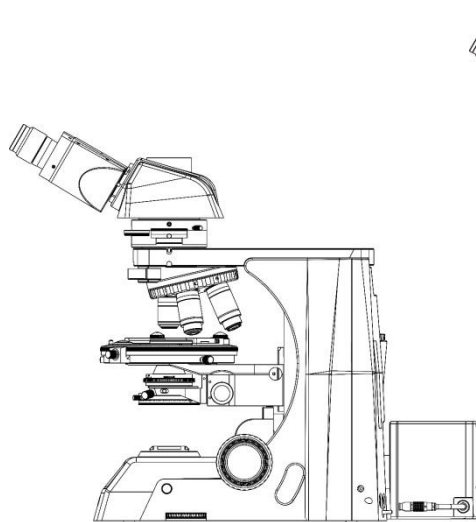


A15.1091-T

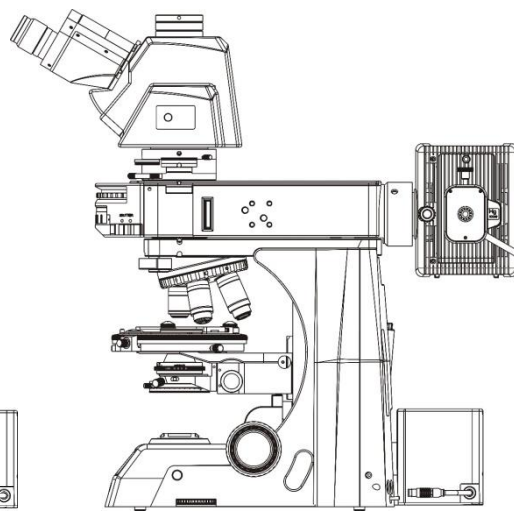

A15.1091-TR

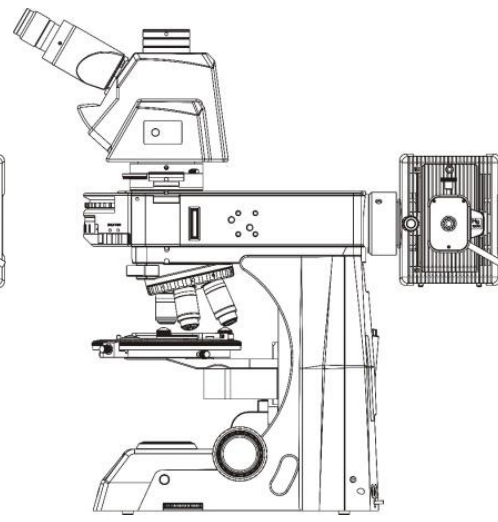

A15.1091-R

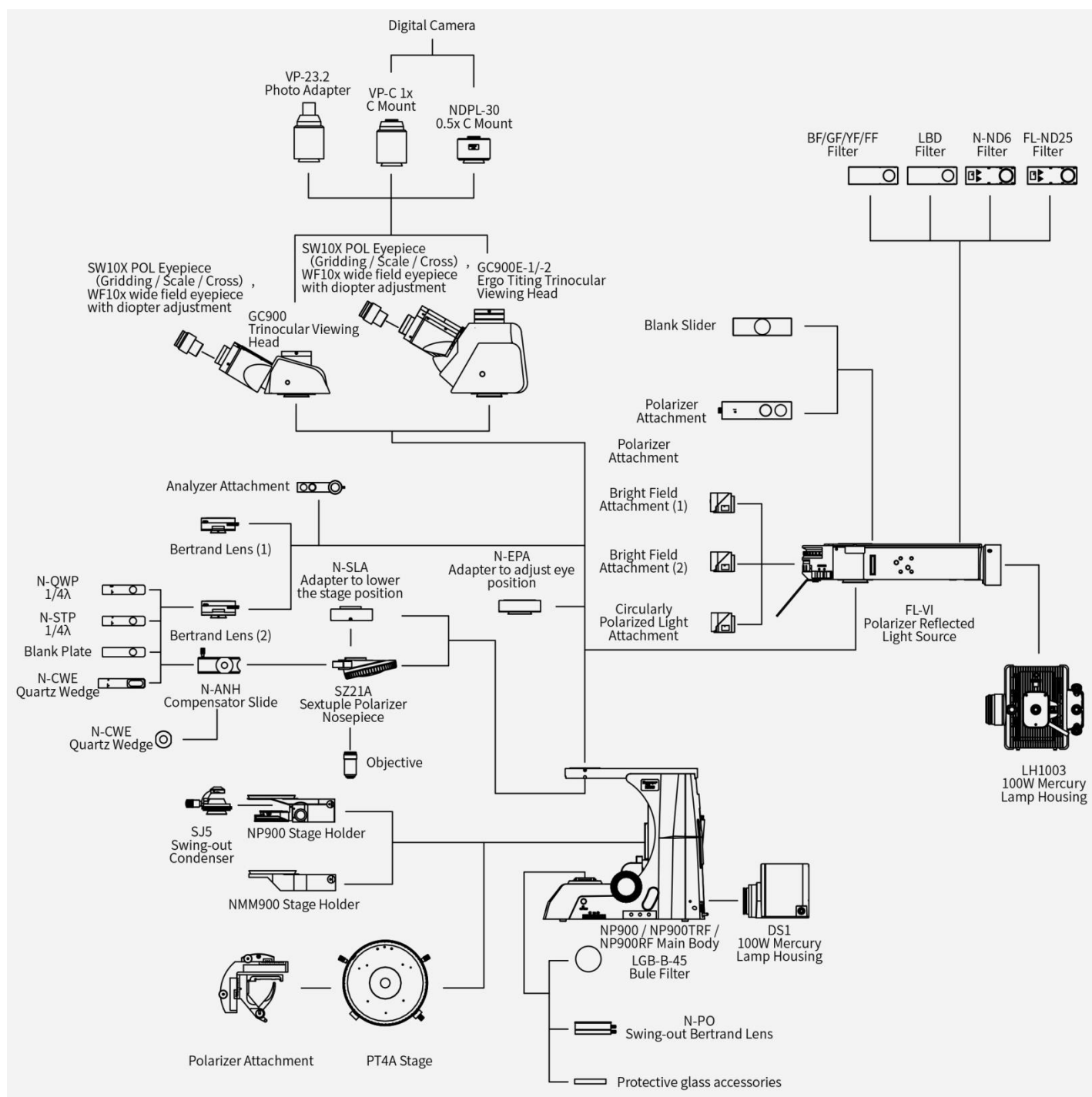

MICROSCOPE IS OUR FOCUS

The Key To Micro World

Opto-Edu is one of the most professional supplier for microscope & educational instruments from China. We have been focusing in this field for more than 25 years, has a professional sales team who can assist our customer to find the BEST SELECTED microscope with competitive price and 3 YEARS Warranty.

25+ Year professional experience we know Microscopes the best!

150+ Microscope & accessory manufacturers supply all models from China

200+ Hot sale microscopes & Newest Models Updated Every Month

750+ Customer from all over the world & keep rising every day

1500+ Microscope products create your one-stop purchase platform

3000+ Educational Instruments For School, College And University Teaching

The Most Professional Microscope Manufacturer in China!



MICRO
WORLD

The Key To OPTO-EDU



Opto-Edu (Beijing) Co., Ltd.

F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China

Official Main Website: www.optoedu.com

www.cnoec.com.cn, www.optoedumicroscope.com, www.microscopemadeinchina.com

Skype: xincnoec Wechat, Mobile, WhatsApp: +86 13911110627

Tel: +86 10 88696085 Email: sale@optoedu.com

