


A14.0960/A16.0960

A14.0960 Research Level Inverted Biological Microscope


A16.0960 Research Level Inverted Fluorescent Microscope

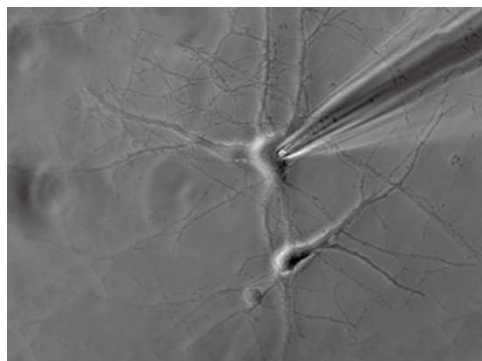

Single-cell Patch Clamp



**Microinjection And Other
Research And Applications**

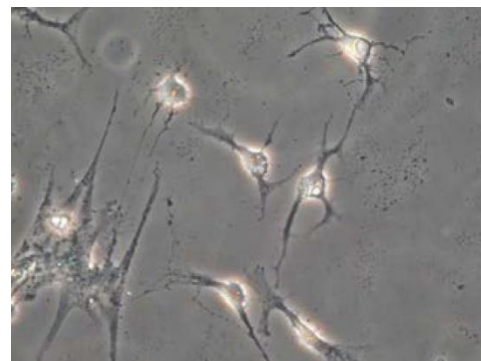

Cell Observation




As a representative of Chinese optical, mechanical, electrical and computing integrated microscope, A16.0960 fully automatic research inverted microscope integrates multiple technologies such as electrical control, fluorescence, phase contrast, relief phase contrast, simple polarization, DIC, and objective lens with correction ring. A three-level electric control of single machine, TPC and computer is realized. A14.0960/A16.0960 has powerful expansion functions, providing a more complete solution for life science research.

With flexible structural design and excellent optical performance, A14.0960/A16.0960 can provide you with high-definition wide-field images, providing powerful technical support for cell observation, single-cell patch clamp, micro-injection and other research and applications.

► Elevation Adjustable Viewing Head

The A14.0960/A16.0960 is equipped with a 20-45° adjustable binocular viewing head, which can raise the eye point by 78mm (65mm interpupillary distance) as required. It can be observed easily and quickly even in a standing state, which is effective avoiding fatigue of eyes and limbs.

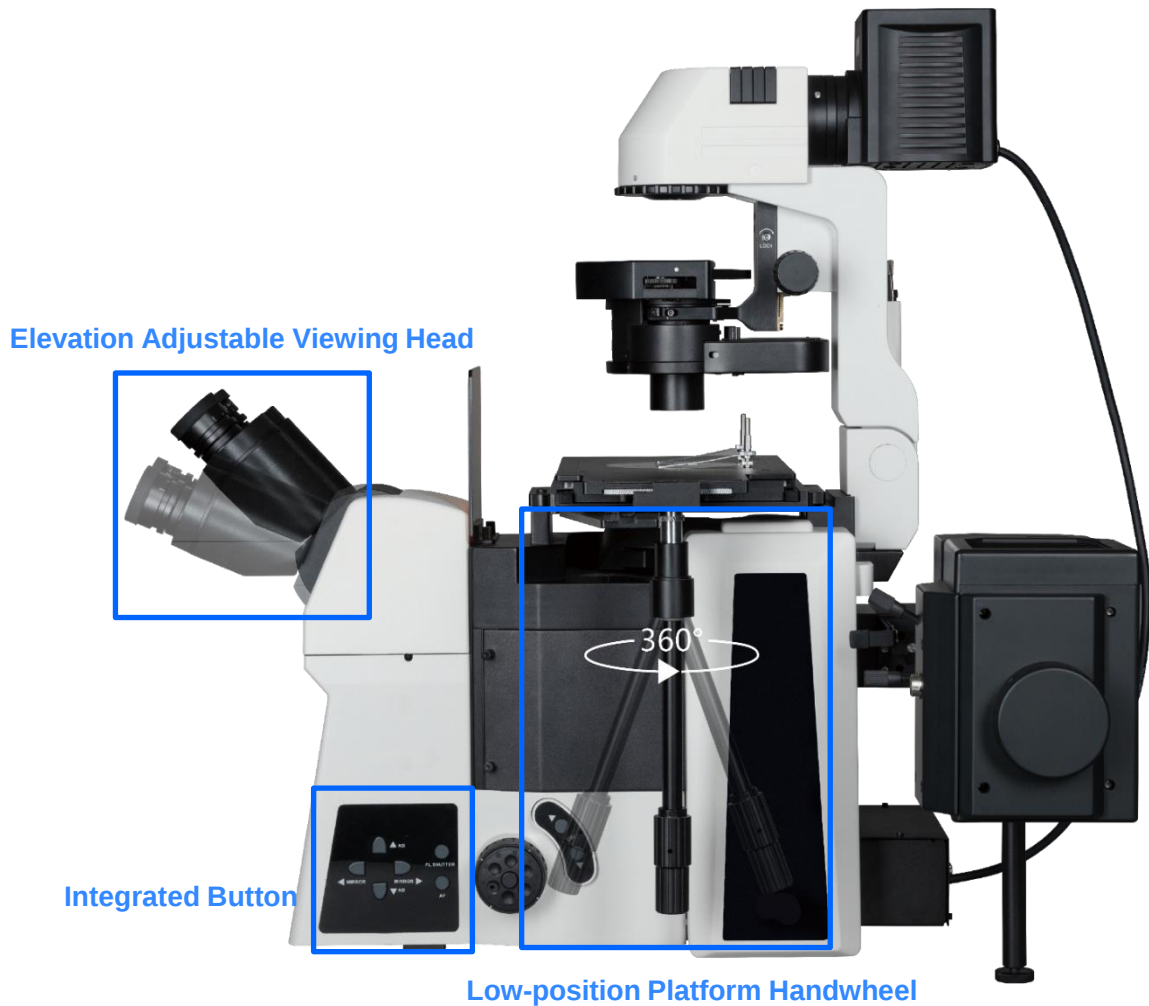

► Integrated Button

A14.0960/A16.0960 retains the traditional coarse and fine adjustment mode, cancels the gear mechanism, integrates electric control technology and realizes manual and automatic integration. The brightness, objectives, attenuator turntable and fluorescent turntable can be quickly switched or rotated through the corresponding buttons on the button panels on both sides. AF one-key auto focus can quickly adjust the Z-axis height according to the real-time image, and eliminates the need for fine-adjusting steps and improves work efficiency.

► Low-position Platform Handwheel

The low-position platform handwheel can be rotated 360°, which can effectively reduce the hand fatigue caused by long-term operation and is ergonomic and improves the convenience of operation.

► Tilting Transmission Lighting Structure

The tilting structure of the transmission bracket ensures a large working space for users and facilitates the replacement of specimens.



► Integrated Panel Control

A14.0960/16.0960 front digital display screen can display the status of multiple components such as microscope objective lens magnification, brightness, and fluorescence band in real time, and can complete the setting of all observation methods, the selection of optical port status, and the switch of fluorescence shutter. Greatly enhance the customer experience and make research work more convenient.



► Waste Liquid Drainage Structure

The drainage groove structure is designed under the objective nosepiece, which can prevent the optical components and modules from accidentally being wetted by the cell culture medium or soaking liquid, contaminating and damaging the microscope, simplifying maintenance.



► High Scalability

The A14.0960/A16.0960 has a very large rack that provides redundant space for third-party configurations. You can choose single, layer or double-layer optical path according to your needs, up to 16 color filters can be replaced; Multi-light port design can be equipped with confocal module, wide, field fluorescent camera, infrared camera, etc., which provides maximum scalability for the research work.


► LWD Semi-APO Series Objectives

Full series of LWD semi-APO objectives were developed for cell culture. The objectives with correction ring have significant advantage for the observation of glass substrates and culture vessels with different thickness, and achieves accurate focusing by correcting the cover glass thickness.


		BF	DF	DIC	POL	FL	PH	RPC	
LWD Infinity Plan Semi-APO Positive Phase Contrast Objective	4X/0.13, WD=17mm,	√			√	√	√		A5C.0960-4
	10X/0.30, WD=8.8mm,	√		√	√	√	√		A5C.0960-10
	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	√	√	√	√	√	√		A5C.0960-20
	40X/0.60, WD=2.85~4.05mm, 0~2mm Cover glass thickness	√	√	√	√	√	√		A5C.0960-40
	60X/0.70, WD=1.42~2.10mm, 0~1.3mm Cover glass thickness	√	√	√	√	√	√		A5C.0960-60
LWD Infinity Plan Semi-APO Objective	4X/0.13, WD=17mm,	√			√	√			A5F.0961-4
	10X/0.30, WD=8.8mm,	√		√	√	√			A5F.0961-10
	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	√	√	√	√	√			A5F.0961-20
	40X/0.60, WD=2.85~4.05mm, 0~2mm Cover glass thickness	√	√	√	√	√			A5F.0961-40
	60X/0.70, WD=1.42~2.10mm, 0~1.3mm Cover glass thickness	√	√	√	√	√			A5F.0961-60
Infinity Plan APOchromatic Objective	4X/0.16, WD=12.8mm	√			√	√			A5F.0962-4
	10X/0.40, WD=3.2mm	√	√	√	√	√			A5F.0962-10
	20X/0.75, WD=0.6mm	√	√	√	√	√			A5F.0962-20
	40X/0.95, WD=0.15mm	√		√	√	√			A5F.0962-40
	60X/0.90, WD=0.26mm	√		√	√	√			A5F.0962-60
	60X/1.25, WD=0.14mm Oil	√			√	√			A5F.0962-60S
	100X/1.35, WD=0.13mm Oil	√			√	√			A5F.0962-100
Infinity Plan Super APO Objective	10X/0.14, WD=3.1mm	√			√	√			A5F.0963-10
	20X/0.80, WD=0.6mm	√			√	√			A5F.0963-20
	40X/0.95, WD=0.18mm	√			√	√			A5F.0963-40
	60X/1.42, WD=0.17mm Oil	√			√	√			A5F.0963-60
	100X/1.45, WD=0.14mm Oil (Under development)	√			√	√			A5F.0963-100
LWD Infinity Plan Semi-APO Relief Hoffman Phase Contrast Objectives	10X/0.30, WD=8.8mm,	√	√	√	√	√		√	A5C.0965-10
	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	√	√	√	√	√		√	A5C.0965-20
	40X/0.45, WD=2.85~4.05mm, 0~2mm Cover glass thickness	√	√	√	√	√		√	A5C.0965-40



Infinity correction plan
APO objectives

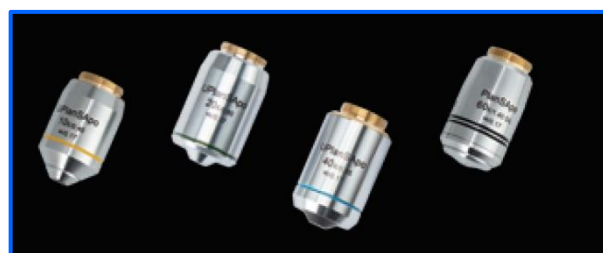

Infinity correction plan
super-APO objectives

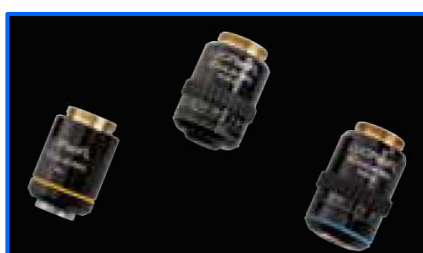

LWD plan semi-APO Relief
phase contrast objectives



LWD plan semi-APO
objectives




LWD plan semi-APO
phase contrast objectives

► Cool LED High-performance Fluorescent Light Source


Compact LED power box, is able to provide illumination to excite common fluorescent filter with different wavelengths, is available to be applied for multi-channel fluorescent filters.




A14.0960/A16.0960 Application Cases

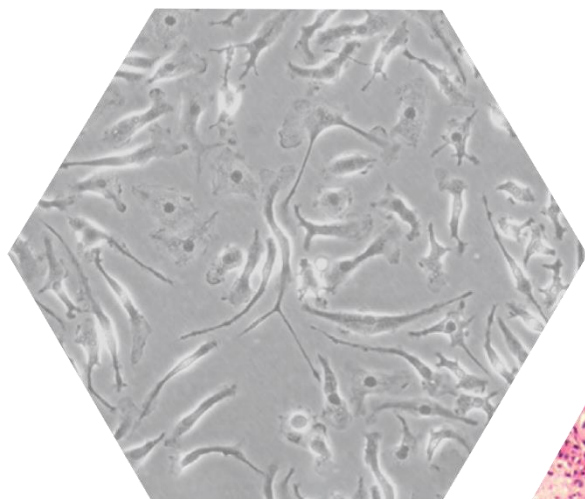


Normal Rat Kidney
Cells Phase Contrast 10X

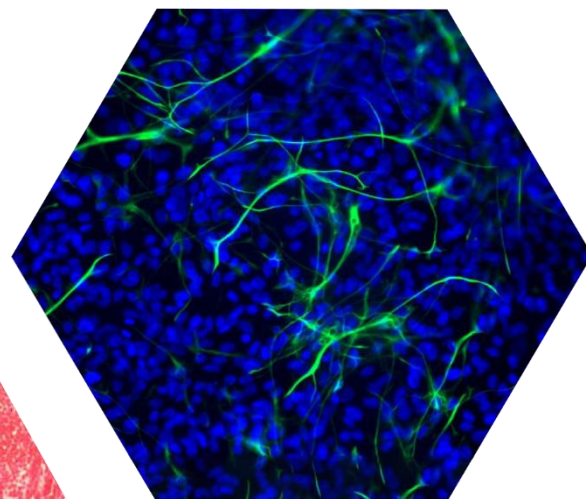


Human Smooth Muscle
Cells Fluorescence 20X

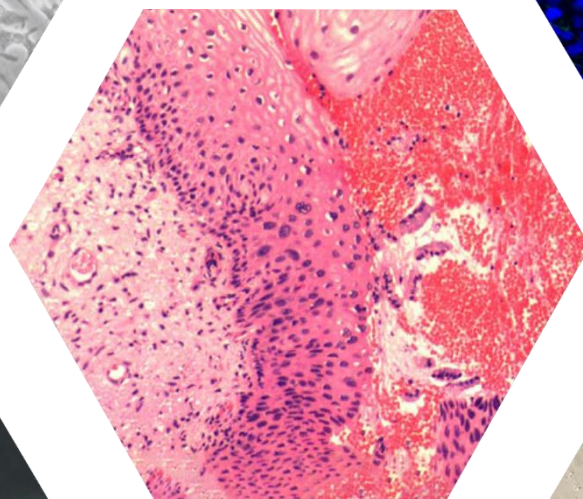

Skin Papilloma
Bright Field 40X

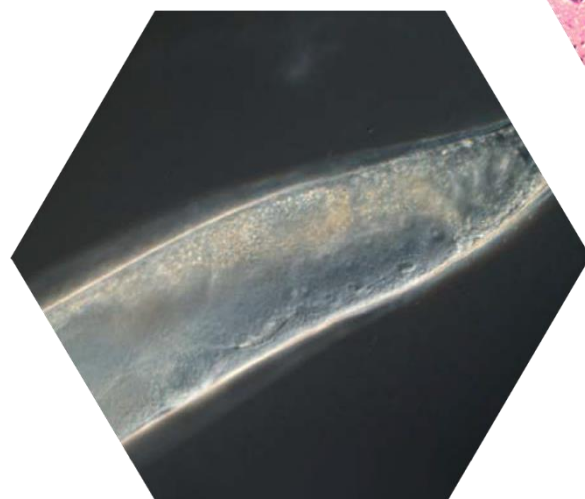

Nematode DIC
40X

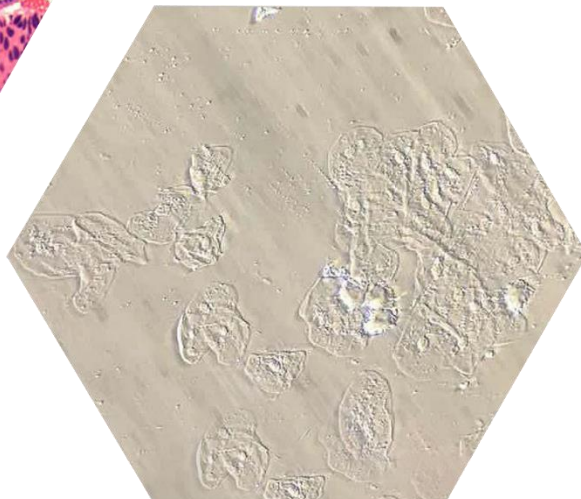

Oral Epithelial Cells
Relief Phase Contrast 10X




A14.0960 Research Level Inverted Biological Microscope A16.0960 Research Level Inverted Fluorescent Microscope			A14.0960	A16.0960
Function	BF, Bright Field		●	●
	PH, Phase Contrast		●	●
	POL, Polarizing		○	○
	DIC, Differential Interfered Contrast		○	○
	Relief Hoffman Phase Contrast		○	○
Head	20°-45° tilting Binocular head, inverted image,interpupillary distance: 50mm~76mm	A53.0960	●	●
Eyeiece	High eyepoint plan eyepiece PL10X/22mm, diopter adjustment	A51.0904-1022T	●●	●●
Nosepiece	Motorized bright filed sextuple nosepiece with DIC slot.			
LWD Infinity Plan	4X/0.13, WD=17mm,	A5C.0960-4	●	●
Semi-APO	10X/0.30, WD=8.8mm,	A5C.0960-10	●	●
Positive Phase	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	A5C.0960-20	●	●
Contrast	40X/0.60, WD=2.85~4.05mm, 0~2mm Cover glass thickness	A5C.0960-40	●	●
Objective	60X/0.70, WD=1.42~2.10mm, 0~1.3mm Cover glass thickness	A5C.0960-60	●	●
LWD Infinity Plan	4X/0.13, WD=17mm,	A5F.0961-4	○	○
Semi-APO	10X/0.30, WD=8.8mm,	A5F.0961-10	○	○
Objective	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	A5F.0961-20	○	○
	40X/0.60, WD=2.85~4.05mm, 0~2mm Cover glass thickness	A5F.0961-40	○	○
	60X/0.70, WD=1.42~2.10mm, 0~1.3mm Cover glass thickness	A5F.0961-60	○	○
Infinity Plan APOchromatic Objective	4X/0.16, WD=12.8mm	A5F.0962-4	○	○
	10X/0.40, WD=3.2mm	A5F.0962-10	○	○
	20X/0.75, WD=0.6mm	A5F.0962-20	○	○
	40X/0.95, WD=0.15mm	A5F.0962-40	○	○
	60X/0.90, WD=0.26mm	A5F.0962-60	○	○
	60X/1.25, WD=0.14mm Oil	A5F.0962-60S	○	○
	100X/1.35, WD=0.13mm Oil	A5F.0962-100	○	○
Infinity Plan Super APO Objective	10X/0.14, WD=3.1mm	A5F.0963-10	○	○
	20X/0.80, WD=0.6mm	A5F.0963-20	○	○
	40X/0.95, WD=0.18mm	A5F.0963-40	○	○
	60X/1.42, WD=0.17mm Oil	A5F.0963-60	○	○
	100X/1.45, WD=0.14mm Oil (Under development)	A5F.0963-100	○	○
Working Stage	Manual mechanical stage, 300mm (X) × 240mm (Y), moving range is 135mm (X) × 85mm (Y),stage thickness is 30mm. Lower universal right side handle. With lock function. During locking, moving range 50mm (X) × 50mm (Y). It can use with 35mm petri dish holder, culture flask holder, sliced holder, Φ110mm replaceable disc stage with clips.	A54.0960	●	●
	Stage: Motorized, 260*153mm with working range 110*75mm, The maximum speed is 3mm/s, the repeat positioning accuracy is ±1um, and it can be adapted to petri dish and slice observation	A54.0963	○	○
Marzhauser Motorized Stage	Marzhauser Motorized stage, moving range: 120mm(X) x 80mm(Y); resolution: 0.02μm, repeatability <1μm, with control box and operating handle, with slider holder.	A54.0961	○	○
	φ36 Petri dish holder, used with Marzhauser-SCAN		○	○
	96-well plate petri dish holder,used with Marzhauser-SCAN		○	○
	Teraseki holder,used with Marzhauser-SCAN		○	○
Prior Motorized Stage	Prior Motorized stage, moving range: 114mm(X) x 75mm(Y), resolution: 0.01μm, repeatability:0.2 μm, with control box and operating handle.	A54.0962	○	○
	φ36 Petri dish holder, used with Prior stage.		○	○
	96-well plate petri dish holder,used with Prior stage.		○	○
	Teraseki holder,used with Prior stage.		○	○
Frame, Focusing & Camera Port	Motorized coarse & fine coaxial focusing, range: 10.5mm, precision: 1 μm.	A54.0960-AV	●	●
	Motorized left camera port, splitting ratio:0:100, 50:50, 100:0.			
	Motorized upper camera port, splitting ratio: 100:0, 0:100,with CTV adapter			
	Dual optical path, with fluorescent light barrier.			
	Base mounting bracket anti-vibration	A54.0960-AV	○	○
	Right camera port, manual two gears switching, splitting ratio: 100:0, 0:100, POV16mm, built-in C-mount 1X. (MUST ordered with whole unit, cannot be added later as optional)	A55.0960	○	○
Remark	●' means standard outfit, 'o' means optional outfit, 'x' means not compatible, 'A','B','C' means only compatible for specific condenser A,B or C			

A14.0960 Research Level Inverted Biological Microscope A16.0960 Research Level Inverted Fluorescent Microscope			A14.0960	A16.0960
Condenser B	Motorized LWD universal 7 positions condenser with N.A.0.55, W.D.27mm. 3 positions ø30mm for Phase contrast observation. 4 positions ø38mm, External connecting cable 100cm Work For BF, PH, DIC.	A56.0962	●	●
Phase Contrast for B	4X phase contrast aperture group	A5C.0962-4	●	●
	10X phase contrast aperture group	A5C.0962-10	●	●
	20X/40X/60X phase contrast aperture group	A5C.0962-204060	●	●
DIC for A,B	DIC slider & Analyzer for transmitted light, Need to select Objectives: A5F.0961-10x,20x,40x,60x A5F.0962-10x,20x,40x,60x A5C.0960-10x,20x,40x,60x	A5C.0964	Bo	Bo
	10X DIC ring, compatible with condenser B	A5C.0964-10	Bo	Bo
	20X DIC ring, compatible with condenser B	A5C.0964-20	Bo	Bo
	40X/60X DIC ring, compatible with condenser B	A5C.0964-4060	Bo	Bo
Condenser C	Manual 5 positions LWD condenser, NA0.30, WD73mm, to meet 4X-60X phase contrast, simple polarization or 10X-40X RPC relief contrast. Work For BF, PH, POL, Relief	A56.0963	o	o
Phase Contrast for C	4X phase contrast aperture group	A5C.0963-4	Co	Co
	10X phase contrast aperture group	A5C.0963-10	Co	Co
	20X/40X/60X phase contrast aperture group	A5C.0963-204060	Co	Co
Polarizing for C	Analyzer, MUST be selected for Reflief Hoffman	A5P.0963-A	Co	Co
	360° rotatable polarizer, MUST be selected for Reflief Hoffman	A5P.0963-P	Co	Co
Relief Hoffman Phase Contrast Slide & Objective	10X aperture group, compatible with condenser C	A5C.0965-10S	Co	Co
	20X aperture group, compatible with condenser C	A5C.0965-20S	Co	Co
	40X aperture group, compatible with condenser C	A5C.0965-40S	Co	Co
	LWD Infinity Plan Semi-APO Relief Hoffman Phase Contrast Objectives	A5C.0965	Co	Co
	10X/0.30, WD=8.8mm,	A5C.0965-10	Co	Co
	20X/0.45, WD=6.5~7.6mm, 0~2mm Cover glass thickness	A5C.0965-20	Co	Co
	40X/0.45, WD=2.85~4.05mm, 0~2mm Cover glass thickness	A5C.0965-40	Co	Co
Transmit Light	Tilting transmitted light system, transmitted Kohler light system, Condenser lifting adjustable, range: 65mm, able to set 4 filters.	A56.0960	●	●
	10W Cool LED Light Box, color temperature 5000K	A56.0960-LED	o	o
	12V100W halogen illumination box, filament center preset.	A56.0960-HAL	●	●
	12V100W halogen lamp PHILIPS 7724	A56.0960-12100	●●	●●
	LBD filter	A56.0960-LBD	●	●
	Green filter	A56.0960-G	●	●
Reflect Fluorescent Mercury Light	L-shaped design reflection fluorescence illuminator.with attenuation plate.	A5F.0960-MER	x	●
	100W Fluorescence mercury light box		x	●
	100W mercury lamp power supply unit (U5 power supply unit)		x	●
	100W Mercury Lamp (OSRAM)		x	●
	5 positions Motorized filter wheel is offered for both excitation and observation, BF/Shutter/ND6/ND25/ND50; External connecting cable 100cm	A5F.0960-FW	x	o
	Extension accessories	A5F.0960-EX	x	●
Fluorescent Filter	Multi-band LED Light Source Module, with four excitation channels,10W high-power LED	A5F.0960-LED4	x	o
	Motorized turret with 8 positions, built-in shutter, External connecting cable 20cm	A5F.0960-60	x	●
	Rotary disk dustproof box	A5F.0960-DB	x	●
	Fluorescence filter(B) (CHROMA)	A5F.0960-B	x	●
	Fluorescence filter(G) (CHROMA)	A5F.0960-G	x	●
	Fluorescence filter(UV) (CHROMA)	A5F.0960-UV	x	●
CAN Card	Fluorescence filter(R1) (CHROMA) EX:AT620/50X, BS: AT655DC, EM: AT690/50M	A5F.0960-R1	x	o
	Need buy computer too to meet electric control, including MV-image-Research	A54.0966-USB	●	●
Computer	Need buy computer too to meet electric control, including MV-image-Research	A54.0966-PCIE	o	o
	Dell Computer Precision T3650(i7 11700/16GB/256GB+1TB/P620) + P2722H	A54.0966-PC	o	o
	7" screen, Motorized control unit	A54.0966-TPC	o	o
Power supply	External connecting cable 200cm, computer connect with power supply box.	A54.0966-EC	●	●
	External system controller.All motorized units are powered by this.		●	●
	Voltage input:90-265V,output:12V100W/24V5A/15V5A/5V5A		●	●
Digital Camera	External connecting cable 200cm, power supply box connect to microscope		●	●
	SONY IMX267 color camera,CCD 8.9M, 32fps@8.9MP,USB3.0,Sensor 1"	A59.0960-8.9M	o	o
	SONY IMX264 color camera,CCD 5.0M, 35fps@5MP,USB3.0,Sensor 2/3"	A59.0960-5.0M	●	●
	SONY ICX695 color camera,CCD 6.0M, 7.5fps@6MP,USB3.0,Sensor 1"	A59.0960-6.0M	o	o
	SONY ICX695 color cold camera,CCD 6.0M, 7.5fps@6MP,USB3.0,Sensor 1"	A59.0960-6.0MC	o	o
	SONY ICX695 mono color camera,6.0M, 7.5fps@6MP,USB3.0,CCD 1"	A59.0960-6.0MM	o	o
	SONY ICX695 mono color cold camera,6.0M, 7.5fps@6MP,USB3.0,CCD 1"	A59.0960-6.0MMC	o	o
Telescope	RISI 9.9M mono color camera(2960x2960), field of view 17.8mm, pixel size 4.25x 4.25um, 30 fps, 73% quantum efficiency, sensor 17.8mm	A59.0960-9.9M	o	o
	Phase Contrast Telescope Φ30	A5C.0910-30	●	●
	Fluorescence Telescope	A5F.0930-2	x	●
C-mount adaptor	C-mount adaptor 0.5X	A55.0930-05	o	o
	C-mount adaptor 0.65X	A55.0930-65	●	●
	C-mount adaptor 1X	A55.0930-10	o	o
Accessories	Allen wrench M3		●●	●●
	Allen wrench M4		●●	●●
	Allen wrench M5		●	●
	Non-fluorescing immersion oil		o	o
Packaging	Carton and foam		●	●
Remark	●' means standard outfit, 'o' means optional outfit, 'x' means not compatible, 'A','B','C' means only compatible for specific condenser A,B or C			


► Microscope revolution

With 50 different image adjustment and processing algorithms, the software is able for photography, image stitching, auto focusing, DOF extension, multichannel acquisition, image measurement. One click to generate a report, reliable and effective.

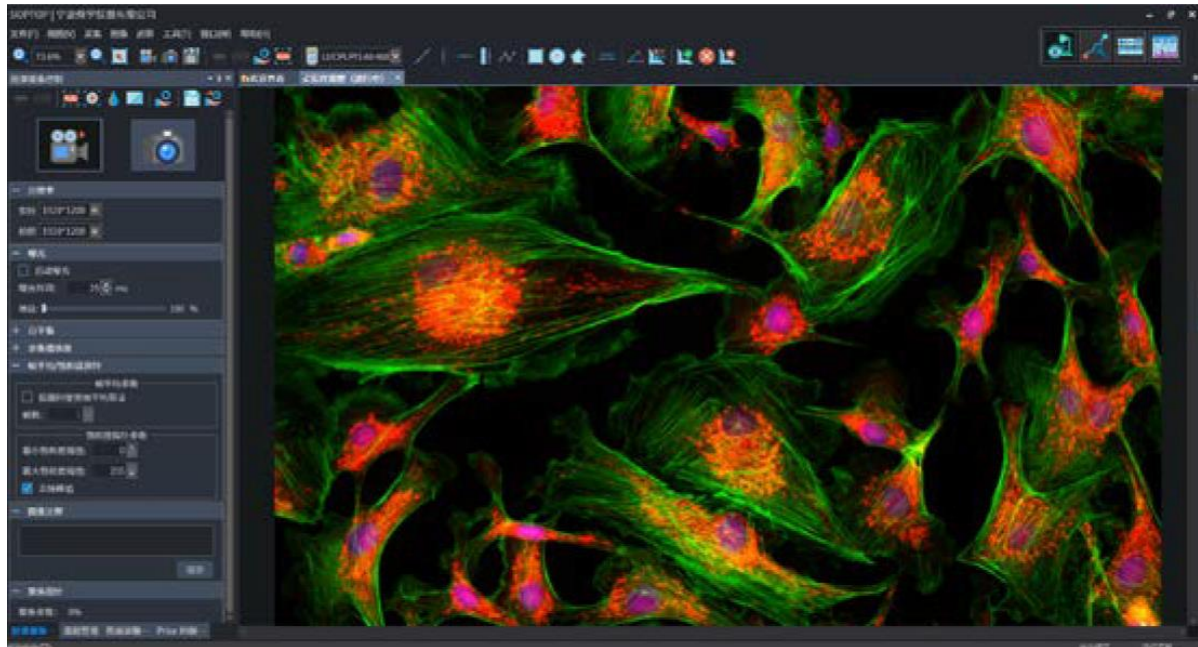


► Image stitching

It is able to stitch the local image, with 8X8 field of view.

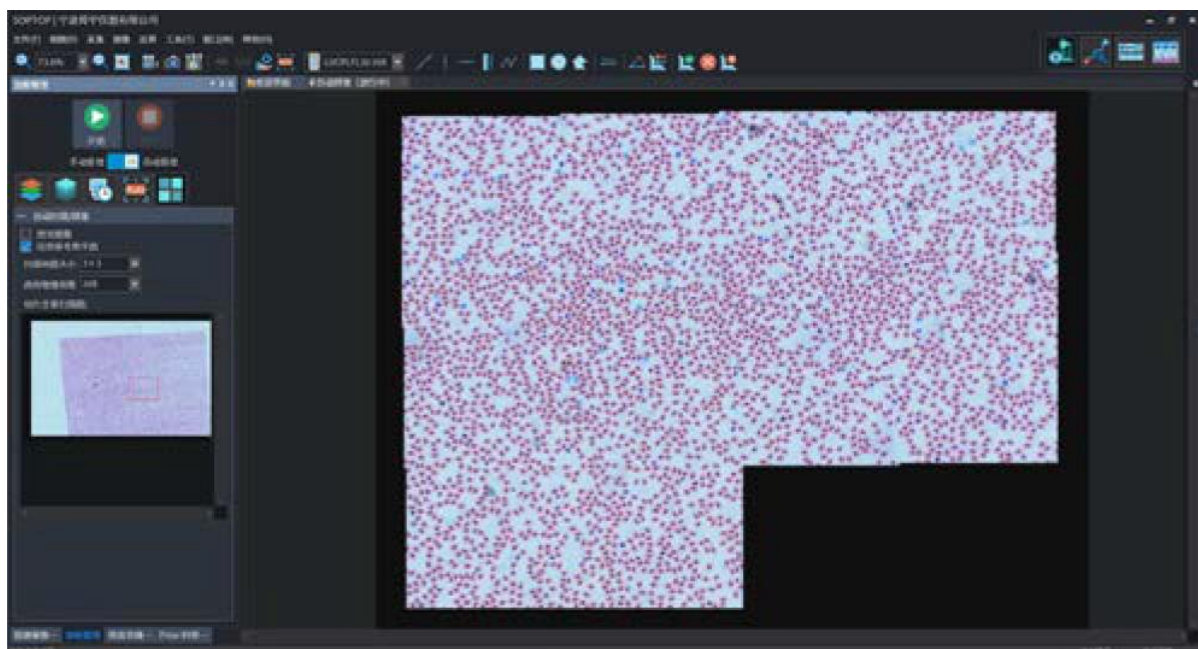


► Cell measurement

Image software with measuring function, is available to analyze cell size, intercellular space, synapse length and so on, including distance, angle, rectangle and round measurement.

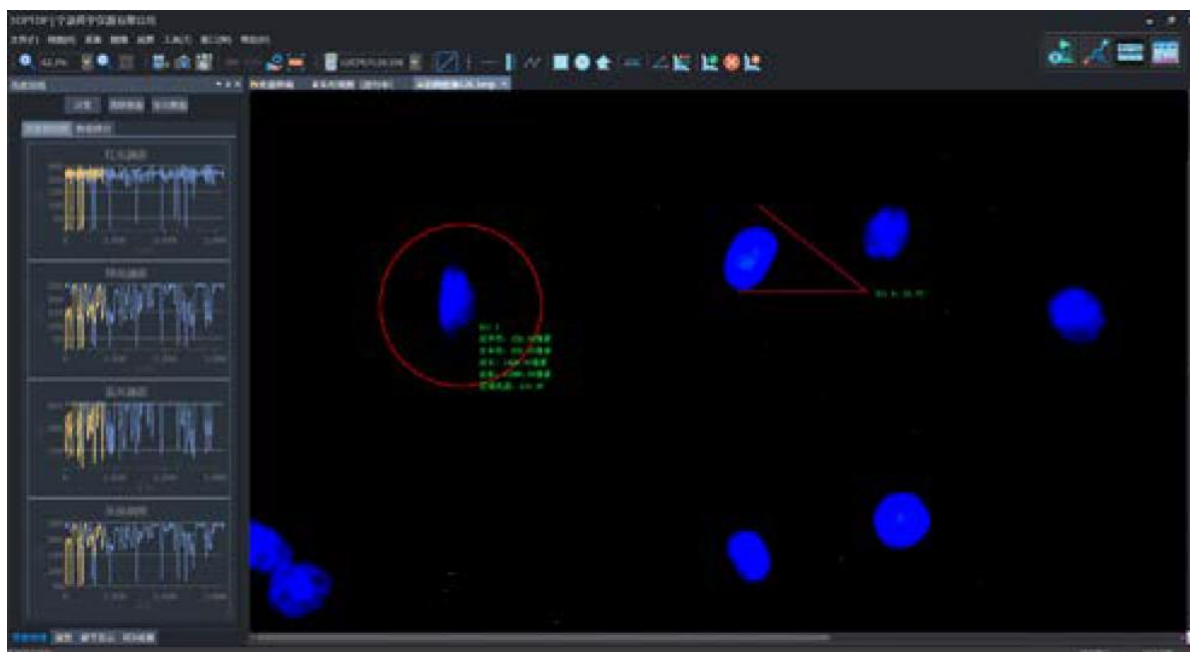


► Fluorescent synthesis

The images from different fluorescent channels, are able to be synthesized automatically and kept separately. Each image is available to be adjusted and the final composited image would be auto updated.

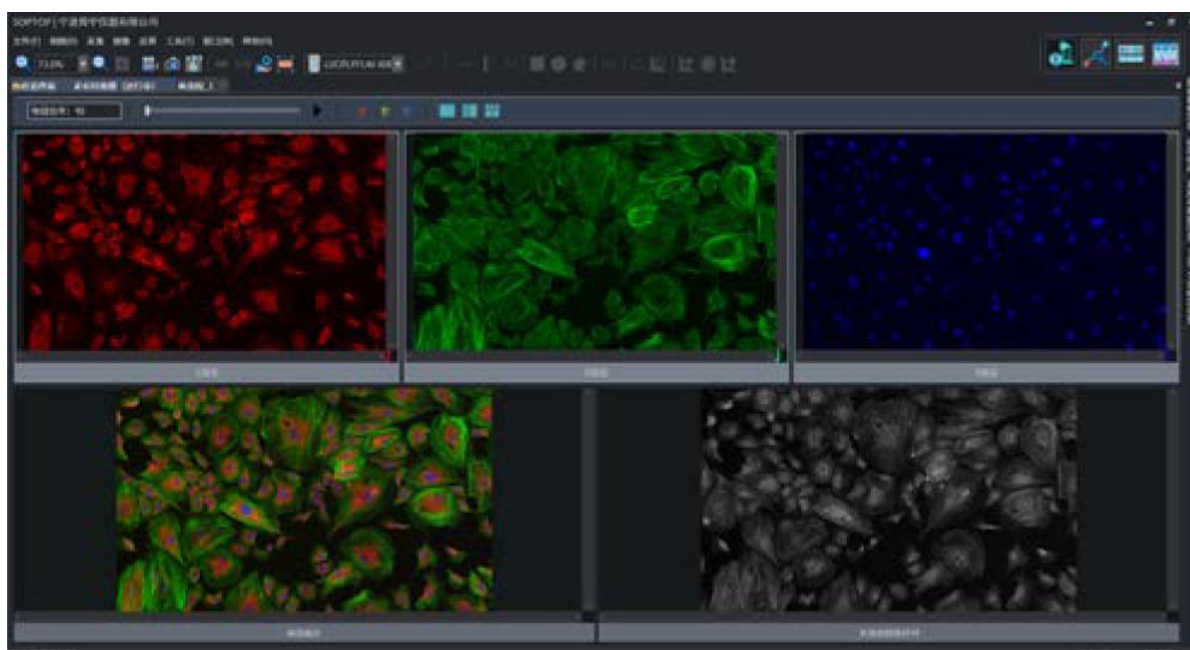


► Report generation

Built-in module to fast generate a report, supports to output WORD document.

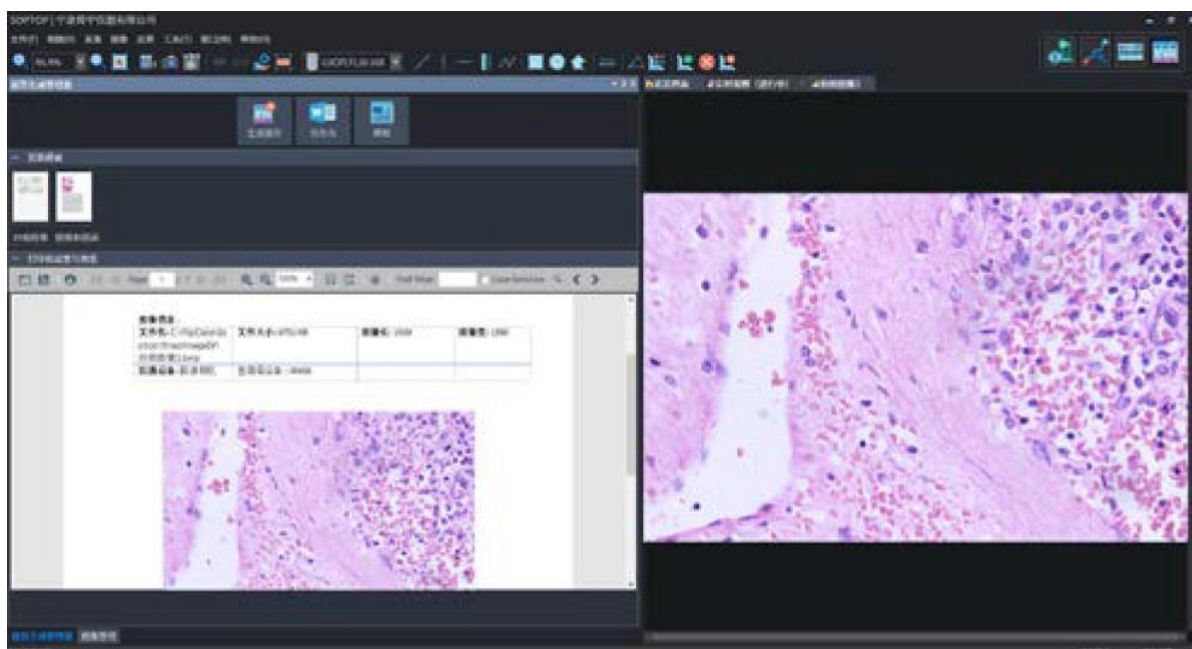

A14.0960/A16.0960 Dimension(mm)

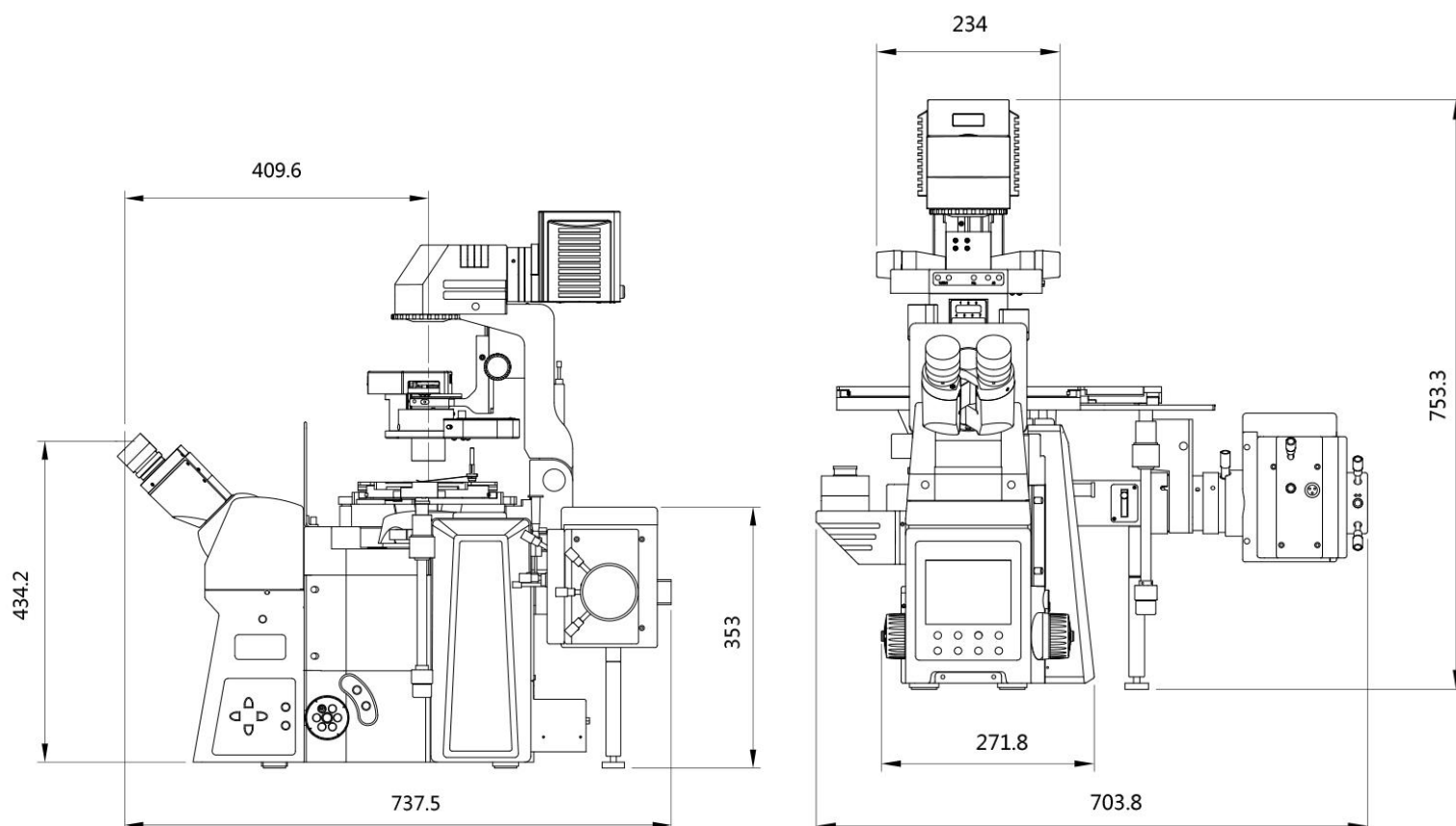



MICROSCOPE IS OUR FOCUS

The Key To Micro World

Opto-Edu is one of the most professional supplier for microscope & educational instruments from China. We have been focusing in this field for more than 25 years, has a professional sales team who can assist our customer to find the BEST SELECTED microscope with competitive price and 3 YEARS Warranty.

25+ Year professional experience we know Microscopes the best!

150+ Microscope & accessory manufacturers supply all models from China

200+ Hot sale microscopes & Newest Models Updated Every Month

750+ Customer from all over the world & keep rising every day

1500+ Microscope products create your one-stop purchase platform

3000+ Educational Instruments For School, College And University Teaching

The Most Professional Microscope Manufacturer in China!



MICRO
WORLD

The Key To OPTO-EDU



Opto-Edu (Beijing) Co., Ltd.

F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China

Official Main Website: www.optoedu.com

www.cnoec.com.cn, www.optoedumicroscope.com, www.microscopemadeinchina.com

Skype: xincnoec Wechat, Mobile, WhatsApp: +86 13911110627

Tel: +86 10 88696085 Email: sale@optoedu.com

